eye
banking
tools

DMEK vacuum donor **punch**
donor **excision** set
slam sutureless amniotic membrane

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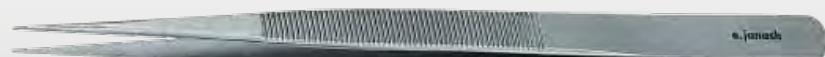
DMEK is the latest accomplishment in deep lamellar surgery. Thanks to the valuable collaboration started with Dr. R. Fogla, our R & D department has been among the first developers of an innovative suction system for this fascinating procedure.

A new project, based on a similar system, developed with the scientific collaboration of the Mestre Eye Bank, has given birth to a new concept of donor punch with peculiar features that facilitate lenticle preparation in Eye Banks.

JD2890FD.100

eyebanking DMEK punch 100 μ calibrated blade

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NEW enhanced donor vacuum trephine. Innovative design allows the surgeon to perform a precise and easy cutting under visual control

DMEK punch con lama $\varnothing$ 10.0 mm, profondità calibrata 100 μ
DMEK punch blade $\varnothing$ 10.0 mm with calibrated depth 100 μ

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DMEK Stripping/peeling forceps for Descemet membrane, delicate anatomically designed tips

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Forceps with tips anatomically designed for the stripping/peeling of the delicate Descemet membrane. The forceps grasp the membrane at the edge and perform the detachment of the Descemet lenticle by means of a peeling action. At this step the membrane is secured to its support by the suction system. Cleavage of the membrane is performed with the Fogla hook J2891E

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DMEK cleavage hook for Descemet membrane

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Delicate blunt hook for Descemet membrane cleavage. The Fogla hook, with atraumatic tips, triangular shaped, performs the cleavage at the edges of the lenticle previously secured to the suction DMEK support. Once the correct thickness of the pure Descemet membrane is achieved, the Janach forceps J2892E performs the detachment by means of the peeling action.

video
www.janach.it

SLAMprocedure sutureless amniotic membrane

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*Blade holder,
designed by M. Campanelli M.D.*

Janach J3713

Lama da Ø 13 mm per patch,
modello M. Campanelli

*Patch blade, Ø 13 mm,
designed by M. Campanelli M.D.*

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*Silicone disc for amniotic membrane,
reusable*

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*Cross-action forceps for holding
and positioning the amniotic
membrane patch*

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Post-operative contact lens applanator

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*Corneal forceps, delicate, 1x2 straight,
0.10 mm teeth, with tying platform, titanium*

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Epithelium removal spatula, for PRK

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Certified
Quality System

UNI EN ISO 9001:2008

UNI CEI EN ISO 13485:2012

© Copyright:
Hermann Janach

e.janach®
www.janach.it

E. Janach S.r.l. Via Borgo Vico, 35 - 22100 Como - Italy
Tel. +39 031 574088 - fax +39 031 572055_international@janach.it

BRO.042.11/14.EN/IT