

DROP^{tonic} 5%

Hypertonic solution for the treatment of corneal edema



- Effective cleansing action of the tear film
- Protective action of the ocular surface
- Greater tolerance



10ml OSD
preservative free
multidose dispenser



Can be used with
contact lenses

Hypertonic Ophthalmics solution

Product Code: DROPTON501VEN; DROPTON501SAG

CE certificate n° QCT- 0079-18 Addendum n° 24-21 Class IIb Notified Body: 0373

CND: Q0299 GMDN: 44237 Volume: 10 ml - 3 ml Shelf Life: 36 months

DESCRIPTION: DROPTonic 5%® is a hypertonic ophthalmic solution for the treatment of corneal oedema occurring as a result of mechanical, toxic, dystrophic or inflammatory damage. By making the tear hypertonic DROPTonic 5%® recalls the excess water from the corneal tissue by osmosis, reducing the oedema-forming state.

COMPOSITION: Riboflavin sodium phosphate (0.000075%), Sodium chloride (5%), Vitamin E TPGS, Sodium Hyaluronate, L-Proline, L-Glycine, L-Leucine, L-Lysine, Dibasic dodecahydrate sodium phosphate, Monohydrate monobasic sodium phosphate, purified water.

CHEMICAL- PHYSICAL PROPERTIES:

- Solution appearance: clear, yellow solution free from particles;
- pH: 6.7-7.3;
- Sterile.

WARNINGS AND PRECAUTIONS: The product is intended for external ophthalmic use only. The product is intended for single patient use. Do not use in cases of known intolerance or hypersensitivity to its components. Do not use in case of damaged or unclosed package before the first opening. Do not use the product after the expiry date indicated on the package. In case of problems during use, discontinue treatment and consult a doctor. During product application, do not touch the eye or any other surface with the dispenser tip. The expiry date refers to unopened and correctly stored product. Before use, carefully read the instructions found on the leaflet. Keep away from children. Do not waste into the environment after use. The solution may stain surfaces and fabrics. In rare cases, corneal micro calcifications have been found in patients with a history of corneal disease treated with ophthalmic preparations containing phosphates. In this case it is recommended to use the product under medical supervision. Transparent straw yellow solution. In case of contact with skin or clothes, is easy to remove. It is possible to report adverse events to the manufacturer via the email address info@iromedgroup.com. By reporting adverse events, you can contribute to the safety of this product. Any severe events should be notified to the manufacturer and the competent authority of the country in which this product has been used.

PRIMARY PACKAGE:

- OSD bottle made for 10 ml in low density polyethylene (LDPE);
- white HDPE eye dropper and a colored cap.


STORAGE: Immediately close the dispenser after use: the remaining product can be used up to 3 months after the first opening. Store away from heat sources and in any case at a temperature less than 25°C.

EXPIRY: With an undamaged package, do not use after the expiry date indicated on both the box and the dispenser. Do not separate the dispenser from the box and the leaflet, to keep all the information for a correct use in the same place.

RECYCLING:

Instructions	Box
PAP 22	PAP 21
PAPER	
Cap	Bottle
LLDPE 4	LDPE 4
PLASTIC	

PRODUCT CODE INFORMATION:

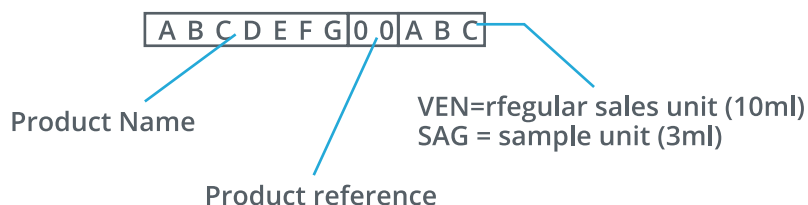

IROMED GROUP S.r.l.
Via Tempio del Cielo 3/5
00144 Rome (Italy)

UK Distributor:
Kestrel Ophthalmics Ltd
01202 696963



SERVImed Industrial S.p.A.
Via Tempio del Cielo 3/5
00144 Rome (Italy)