




The future of vision

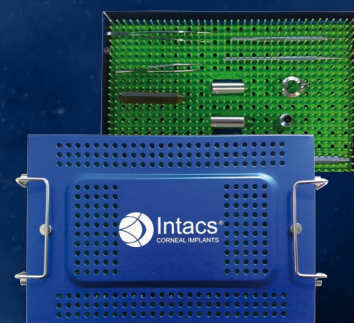

Intacs® & Intacs® SK Corneal Implantation Instruments

For use with both Intacs® & Intacs® SK Corneal Implants				
<i>Product Description</i>	Intacs® Surgical Instruments consist of a series of ophthalmic instruments that are designed for the placement of Intacs® Corneal Implants.			
<i>Indications</i>	Intacs® Instruments are indicated for use in patients undergoing ophthalmic surgery requiring a lamellar dissection of the cornea. They can be used for Intacs® (optical zone 7mm) and Intacs® SK (optical zone 6mm).			
<i>Presentation</i>	Supplied non-sterile.			
				
<i>Reference</i>	KV10445	KV10478-2	KV10606	KV1160
<i>Instrument</i>	Intacs® Forceps	Vacuum Centering Guide, Sloped	Pocketing Hook	Symmetric Glide
<i>Material</i>	Titanium	Stainless steel	Titanium	Titanium
<i>Characteristics</i>	- Notched tip. - Facilitates easy removal of the segments from the carrier.	- Port for attachment to a vacuum source. - Equipped with fins on the bottom. - Centers and stabilizes the Corneal Separators onto the cornea during the dissection of the intrastromal tunnels.	- Alloy probe. - Initiates the delamination of the stroma.	- Bifurcated glide. - Guides in both the CW and CCW Corneal Separators.

CE

REV: 03

For use only with Intacs® Corneal implants (optical zone 7mm)



Reference	KV10730-1	KV10730-2	KV10732	KV10733
Instrument	Corneal Separator, CW, 0.9 mm		Inspection Gauge, 0.9 mm	Procedure Marker, 0.9 mm
Material	Stainless steel		Ultem	Stainless steel
Characteristics	· Also known as “Corneal Dissectors”. · 210° blade attached to the body by a spoke. · Creates a CW and CCW annular subsurface tunnel.		· Also known as “Separator or Dissector Inspection Gauge”. · Inspects the concentricity and planarity of the 0.9mm Corneal Separator blades.	· Also known as “Incision & Placement Marker”. · Marks the incision. · Aligns the Vacuum Centering Guide.

For use only with Intacs® SK Corneal Implants (Optical Zone 6mm)



Reference	AT11209-1	AT11209-2	AT11210	AT11211
Instrument	SK Corneal Separator, CW 0.6 mm		SK Inspection Gauge	Procedure Marker
Material	Stainless steel		Ultem	Stainless steel
Characteristics	· Also known as Corneal Dissectors. · 210° blade attached to the body by a spoke. · Creates a CW and CCW annular subsurface tunnel.		· Also known as “Separator or Dissector Inspection Gauge”. · Inspects the concentricity and planarity of the 0.6mm Corneal Separator blades.	· Also known as “Incision & Placement Marker”. · Marks the incision. · Aligns the Vacuum Centering Guide.

Additional Instruments

Description and References

- Diamond Knife, 15° angled blade - Ref: ME410-P01.
- Sterilization Tray - Ref: 4710-A.
- Sinskey Hook - Ref: 6-251.
- Bonn Suturing Forceps - Ref: 2-110.
- Ring Marker with Cross Hairs, 11mm Zone - Ref: 9-789W-1.
- Kratz Barraquer Speculum - Ref: 9-556F.
- Vacuum system - Ref: KV2000.
- Vacuum tubings - Ref: 7K5000.