



INTERNATIONAL

FORM FIT®

INTRACANALICULAR PLUGS

LONG TERM OCCLUSION


FORM FIT® is truly a one size plug solution. Once inserted into the canaliculus, FORM FIT® hydrates in approximately 10 minutes. As it hydrates, FORM FIT® increases in size until it completely fills the vertical canalicular cavity.

HOW IT WORKS




FORM FIT® is made of a safe and biocompatible hydrogel material that expands into a soft, pliable, gelatinous material when it contacts tear film. FORM FIT® securely conforms to the vertical canaliculus so there is no need for a cap on the end of the plug. This greatly increases patient comfort and decreases foreign body sensation.

PRODUCT FEATURES

- Hydrates in approximately 10 minutes
- Easily removed via irrigation of the tear duct
- Biocompatible, no foreign body sensation
- Sterile
- Conveniently pre-loaded in inserter
- 0.3mm, one size fits all

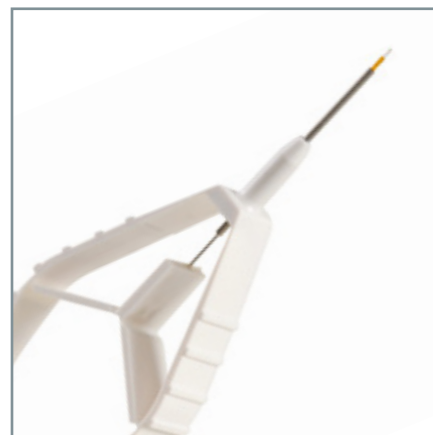

*CE Marked

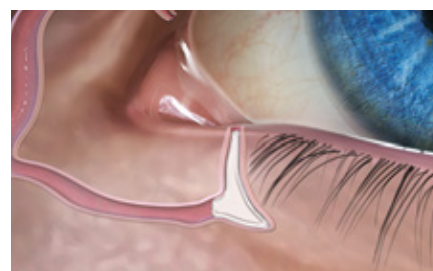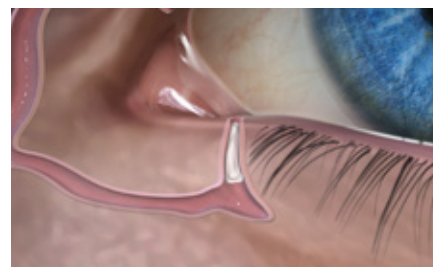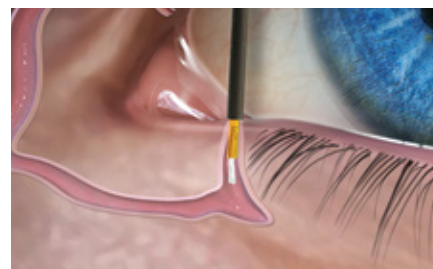

Product Name	Description	Diameter	Ref #
FORM FIT Intracanalicular Plug	10pk	0.3mm	6303
FORM FIT Intracanalicular Plug	2pk	0.3mm	6303-T






SILICONE PUNCTUM PLUGS

LONG TERM OCCLUSION

The distinctive low-profile cap combined with a softer, flexible silicone maximizes patient comfort. Now available with two inserter options. SOFT PLUG® Silicone Punctum Plugs are removable with forceps per instruction for use.

PRODUCT FEATURES

- **NEW** Pen-Style Inserter with dilator tip - 6pks
- Original Inserter with dilator tip - 2pks, 6pks, 20pks
- Sterile, preloaded on inserter in 2-pack or 6-pack dispenser box
- Low-profile cap intended to reduce dislocation
- Economy pack includes 2 inserters and 20 plugs. Holders equipped with an alignment slot for quick loading

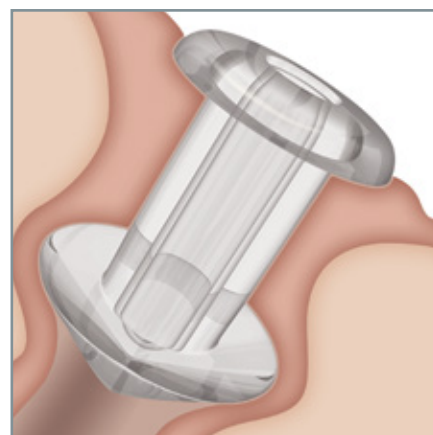

*CE Marked



Pen-Style Inserter
*Pending CE Mark



Original Inserter

Product Name	Description	Diameter	Ref #
SOFT PLUG Silicone Punctum Plug	Micro - 2pk	0.4mm	6610
SOFT PLUG Silicone Punctum Plug	Mini - 2pk	0.5mm	6611
SOFT PLUG Silicone Punctum Plug	Petite - 2pk	0.6mm	6612
SOFT PLUG Silicone Punctum Plug	Small - 2pk	0.7mm	6613
SOFT PLUG Silicone Punctum Plug	Medium - 2pk	0.8mm	6614

Product Name	Description	Diameter	Original Inserter Ref #	Pen-Style Inserter Ref #
SOFT PLUG Silicone Punctum Plug	Micro - 6pk	0.4mm	6610-D	NEW! 7710-D
SOFT PLUG Silicone Punctum Plug	Mini - 6pk	0.5mm	6611-D	NEW! 7711-D
SOFT PLUG Silicone Punctum Plug	Petite - 6pk	0.6mm	6612-D	NEW! 7712-D
SOFT PLUG Silicone Punctum Plug	Small - 6pk	0.7mm	6613-D	NEW! 7713-D
SOFT PLUG Silicone Punctum Plug	Medium - 6pk	0.8mm	6614-D	NEW! 7714-D

Product Name	Description	Diameter	Ref #
SOFT PLUG Silicone Punctum Plug	Micro - 20pk	0.4mm	6610-EP
SOFT PLUG Silicone Punctum Plug	Mini - 20pk	0.5mm	6611-EP
SOFT PLUG Silicone Punctum Plug	Petite - 20pk	0.6mm	6612-EP
SOFT PLUG Silicone Punctum Plug	Small - 20pk	0.7mm	6613-EP
SOFT PLUG Silicone Punctum Plug	Medium - 20pk	0.8mm	6614-EP






SILICONE PUNCTUM PLUG, FLOW CONTROL

LONG TERM OCCLUSION WITH PARTIAL DRAINAGE

The distinctive low-profile cap combined with a softer, flexible silicone maximizes patient comfort. The inserter minimizes potential for accidental release of the plug. Plug is designed with a center channel allowing for partial tear flow. SOFT PLUG® Flow Control Plugs are removable with forceps (per instruction for use).

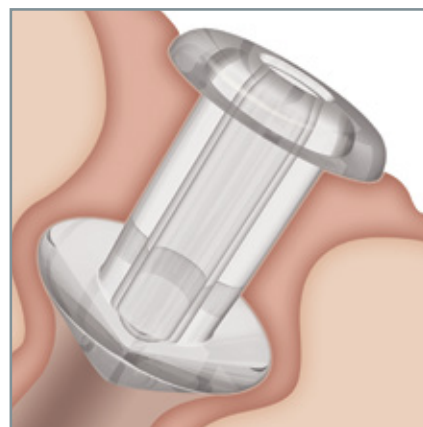

*CE Marked

PRODUCT FEATURES

- 2 plugs per box
- Silicone
- Sterile, preloaded on inserter
- For partial occlusion with limited tear drainage


Product Name	Description	Diameter	Ref #
SOFT PLUG Silicone Punctum Plug, Flow Control	Petite - 2pk	0.6mm	6612-FC
SOFT PLUG Silicone Punctum Plug, Flow Control	Small - 2pk	0.7mm	6613-FC
SOFT PLUG Silicone Punctum Plug, Flow Control	Medium - 2pk	0.8mm	6614-FC



COLLAGEN INTRACANALICULAR PLUGS

TEMPORARY OCCLUSION.

The Collagen Plugs have a smooth proprietary finish that simplifies insertion of the plug. Collagen Plugs are ideal for short-term, diagnostic, and postsurgical occlusion.

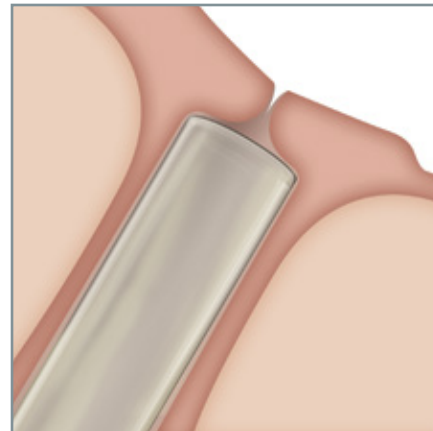

PRODUCT FEATURES

- Absorbable material lasts 2-5 days
- Proprietary finish simplifies insertion
- Evaluate patient acceptance
- 2mm length
- Sterile
- Insert with forceps

Product Name	Qty	Diameter	Ref #
SOFT PLUG Collagen Plug	60 per box	0.2mm	6502
SOFT PLUG Collagen Plug	60 per box	0.3mm	6503
SOFT PLUG Collagen Plug	60 per box	0.4mm	6504



MEDIUM TERM OCCLUSION FOR UP TO 3 OR 6 MONTHS

SOFT PLUG® Extended Duration Punctum Plugs are absorbable canalicular occluders that are an effective option for transient dry eye symptoms before or after surgery, seasonally, etc. and for patients seeking to experience non-permanent options of punctal occlusion.


PRODUCT FEATURES

- Pre and post-surgical use
- Sterile
- Insert with forceps
- 2mm length

SOFT PLUG® EXTENDED DURATION

180 CANALICULAR PLUGS

Up to 180 days.

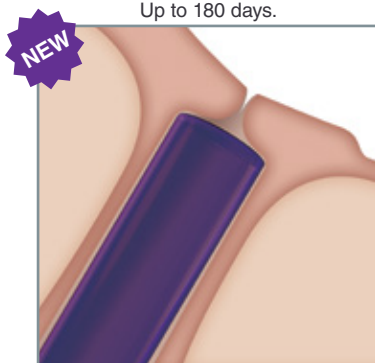


*CE Marked

SOFT PLUG® EXTENDED

DURATION PLUGS

Up to 90 days.

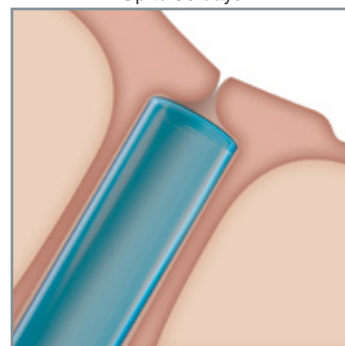


*CE Marked

Product Name	Description	Diameter	Ref #
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug	20pk	0.3mm	6203-A
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug	20pk	0.4mm	6204-A
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug	20pk	0.5mm	6205-A
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug	40pk	0.3mm	6203
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug	40pk	0.4mm	6204
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug	40pk	0.5mm	6205
NEW! SOFT PLUG Extended Duration 180 Canalicular Plug, Variety Pack		3 pair 0.3mm, 4 pair 0.4mm, 3 pair 0.5mm	6200V
SOFT PLUG Extended Duration Plug	20pk	0.2mm	6402-A
SOFT PLUG Extended Duration Plug	20pk	0.3mm	6403-A
SOFT PLUG Extended Duration Plug	20pk	0.4mm	6404-A
SOFT PLUG Extended Duration Plug	20pk	0.5mm	6405-A
SOFT PLUG Extended Duration Plug	40pk	0.2mm	6402
SOFT PLUG Extended Duration Plug	40pk	0.3mm	6403
SOFT PLUG Extended Duration Plug	40pk	0.4mm	6404
SOFT PLUG Extended Duration Plug	40pk	0.5mm	6405
SOFT PLUG Extended Duration Plug, Variety Pack		1 pair 0.2mm, 3 pair 0.3mm, 4 pair 0.4mm, 2 pair 0.5mm	6400V

ACCESSORIES & DRY EYE DIAGNOSTICS



***CE Marked**

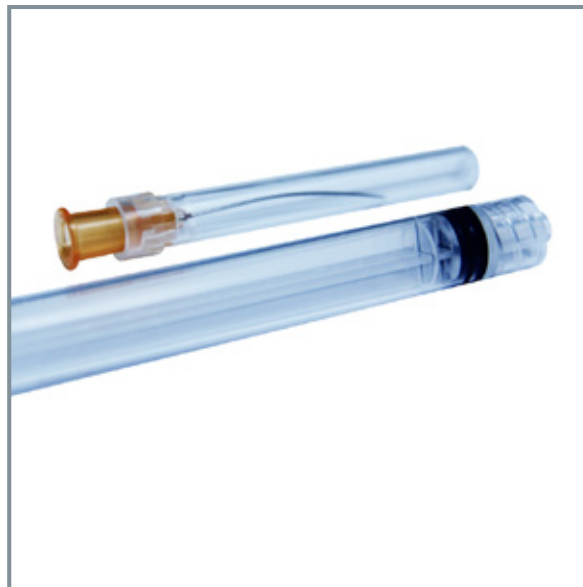

***CE Marked**

Product Name	Description	Qty	Ref #
Punctal Sizing Gauge*	0.5mm/0.6mm	1 per box	6625
Punctal Sizing Gauge*	0.7mm/0.8mm	1 per box	6627
Punctal Sizing Gauge Set*	Includes both instruments and sterilization tray	1 set per box	6630

***CE Marked**

Product Name	Description	Qty	Ref #
Straight Jeweler's Forceps*	Straight #3	1 per box	6643
Curved Jeweler's Forceps*	Curved #7	1 per box	6647
Grooved Forceps*	For the insertion of Extended Duration & Collagen Plugs	1 per box	6648
Castroviejo Lacrimal Dilator*	Double Ended	1 per box	6640

***CE Marked**

Product Name	Description	Qty	Ref #
Lacrimal Irrigation Kit*	One 3ml syringe and one 25g x 1 1/8" lacrimal cannula, curved, open end	10 per box	4003-IK

***CE Marked**





HOT & COLD EYE MASK

RELIEF FOR DRY EYES

Oasis REST & RELIEF® provides heat therapy for the relief of dry eye symptoms. Use as directed.



*Pending CE Mark

PRODUCT FEATURES

- Temperature retaining beads that conform to the body in a soft knit cloth mask for maximum comfort
- Adjustable strap for comfortable compression

Product Name	Description	Qty	Ref #
Oasis REST & RELIEF EYE MASK	Eye mask for hot and cold use	1 mask per box	OM8000

MARKETING MATERIALS

OASIS® offers educational materials in the form of brochures, displays, posters, and videos for use in your practice and to distribute to your patients.

To view and download digital versions of all resources, go to Practice & Patient Materials at www.oasismedical.com

	Product Name	Ref #
	Surgical Catalog	LIT-SRG-CAT
	Dry Eye Symptoms Wall Chart	LIT-DES-WC

	Product Name	Ref #
	Punctal Plug Wall Chart	LIT-WCDE
	OASIS Iris Expander Wall Chart	LIT-WCIE1013



"I usually place the SOFT PLUG Silicone Punctal Plugs. I tell patients that this is a long-term solution that can be reversed. They like hearing that, because most often these patients are younger and are hesitant to do something that may be permanent."

Francis A. D'Ambrosio, Jr, MD.



Glendora, CA, USA

(844) 820-8940 (USA Toll Free) | (909) 305-5400

international@oasismedical.com | www.oasismedical.com

LIT-DEInt-CAT Rev0 12/2018